KMINA

USER MANUAL
MANUAL DE USUARIO
BENUTZERHANDBUCH
MANUEL D'UTILISATION
MANUALE UTENTE
ユーザーマニュアル

K10094 K10095

KMINA - LIGHTWEIGHT SCOOTER
KMINA - SCOOTER LIGERO
KMINA - LEICHTER SCOOTER
KMINA - SCOOTER LÉGER
KMINA - SCOOTER LEGGERO
KMINA - LICHTE SCOOTER
KMINA - 全地形対応スクーター

CE

1 COMPONENTS / COMPONENTES / KOMPONENTEN / COMPOSANTS / COMPONENTI / ONDERDELEN / コンポーネント

	ENGLISH	ESPAÑOL	DEUTSCH	FRANÇAIS	ITALIANO	NEDERLANDS	日本語
1.	Handlebar	Manillar	Lenker	Guidon	Manubrio	Stuur	ハンドルバー
2.	Control panel	Panel de control	Bedienfeld	Panneau de commande	Pannello di controllo	Bedieningspaneel	操作パネル
3.	Lock lever	Palanca bloqueo manillar	Verriegelungs-hebel	Levier de verrouillage	Leva di bloccaggio	Vergrendelings-hendel	ロックレバー
4.	Throttle lever	Palanca de aceleración	Gashebel	Levier d'accélération	Levier d'accélération	Gashendel	アクセルレバー
5.	Front light	Luz delantera	Frontlicht	Feu avant	Luce anteriore	Voorlicht	フロントライト
6.	Front wheels	Ruedas delanteras	Frontlicht	Roues avant	Ruote anteriori	Voorwielen	前照灯
7.	Footrests	Reposapiés	Fußstütze	Repose-pieds	Poggiapiedi	Voetsteun	フットレスト
8.	Rear wheels	Ruedas traseras	Hinterräder	Roues arrière	Ruote posteriori	Achterwielen	後輪
9.	Anti-tip wheels	Ruedas antivuelco	Kippschutzräder	Roues anti-basculement	Ruote antiribaltamento	Anti-kantelwielen	転倒防止輪
10.	Folding supports	Soportes plegado	Klappstützen	Supports de pliage	Supporti di piegatura	Inklapsteunen	折りたたみサポート
11.	Removable battery	Batería extraíble	Herausnehmbarer Akku	Batterie amovible	Batteria estraibile	Uitneembare accu	取り外し可能バッテリー
12.	Seat	Asiento	Sitz	Siège	Sedile	Zitting	座面
13.	Cushion	Cojín	Kissen	Coussin	Cuscino	Kussen	クッション
14.	Frame	Estructura	Rahmen	Châssis	Telaio	Frame	フレーム
15.	Seat folding lever	Palanca plegado del asiento	Sitz-Klapphebel	Levier de pliage du siège	Leva di chiusura del sedile	Hendel voor het inklappen van de zitting	座面折りたたみレバ
16.	Flip-up armrest	Reposabrazos abatible	Hochklappbare Armlehne	Accoudoir rabattable	Bracciolo pieghevole	Opklapbare armleuning	跳ね上げ式アームレスト
17.	Fold-down backrest	Respaldo abatible	Umklappbare Rückenlehne	Dossier rabattable	Schienale ribaltabile	Neerklapbare rugleuning	折りたたみ式背もたれ
18.	NFC card (x2)	Tarjeta NFC (x2)	NFC-Karte (x2)	Carte NFC (x2)	Scheda NFC (x2)	NFC-kaart (x2)	NFCカード (x2)
19.	Brake lever	Palanca de freno	Bremshebel	Levier de frein	Leva del freno	Remhendel	ブレーキレバー
20.	Charger	Cargador	Ladegerät	Chargeur	Caricabatterie	Oplader	充電器

2 MEASURES / MEDIDAS / MESSUNGEN / MESURES / MISURE / AFMETINGEN / 測定方法

WIDTH	WIDTH WITHOUT ARMRESTS	HEIGHT	LENGTH	FOLDING DIMENSIONS
55 cm (21.7”)	44 cm (17.3”)	95 cm (37.5”)	99 cm (39”)	99x44x40 cm (39”x17”x15.7”)
				

SEAT HEIGHT	SEAT DIMENSION	BACKREST HEIGHT	Ø REAR WHEEL	Ø FRONT WHEEL
52 cm (20.5”)	42 x 35 cm (16.5”x14”)	39,5 cm (15.5”)	20 x 4 cm (8”x1.6”)	18 x 3,5 cm (7”x1.4”)

MAX. USER WEIGHT	WEIGHT WITH BATTERY	WEIGHT WITHOUT BATTERY	WEIGHT WITHOUT BATTERY AND ARMRESTS
120 kg (265 lb)	17,7 kg (39 lb)	15,8 kg (35 lb)	15 kg (33 lb)

EN | The KMINA - Lightweight scooter has been designed to offer comfort and safety to the user when moving around. It supports a maximum weight of 265 lb and is designed for independent patients, elderly people or people with reduced mobility.

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3 ASSEMBLY / MONTAJE / VERSAMMLUNG / MONTAGE / MONTAGGIO / MONTAGE / 組立

A - HANDLEBAR / MANILLAR / LENKERWINKELS / POIGNÉE / MANUBRIO / STUUR / ハンドルバー

EN | Loosen the lock lever on the handlebar and raise the handlebar to the most comfortable angle for the user. Tighten the lock lever again until the handlebar is fully secured.

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B - ASIENTO / SEAT / SITZ / SIÈGE / SEDILE / ZITTING / シート

EN | Press the seat folding lever and lift the seat until the frame is fully upright. Fold the fold-down backrest fully upward until it is completely extended, then place the cushio**n.**

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EN | Insert the armrests into the side tubes located behind the backrest. To remove them, press the button on the side tubes and pull the armrests outward.

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EN | 1. To fold the scooter, press the seat folding lever located under the seat and slowly push the backrest and seat downward.
2. Loosen the lock lever on the handlebar and slowly push the handlebar downward.
3. If you want to store the folded scooter upright, the folding supports allow it to stand on its own.
4. To unfold the scooter, repeat steps 1 and 2 in reverse order.

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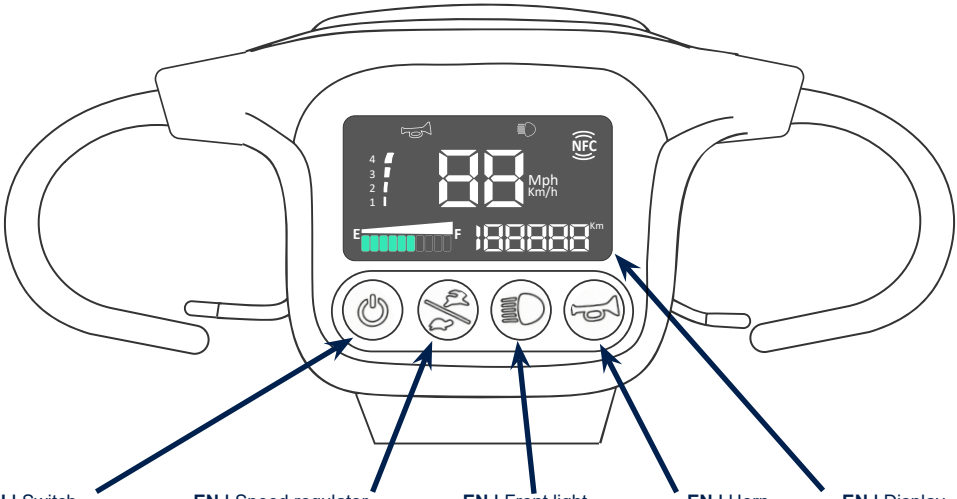
*

EN |

- To switch between metric and imperial units: with the scooter unlocked, press and hold the speed control and the front light for 5 seconds until you hear a beep.
- To enable/disable the reverse warning sound: with the scooter unlocked, press and hold the front light and the horn for 5 seconds until you hear a beep.

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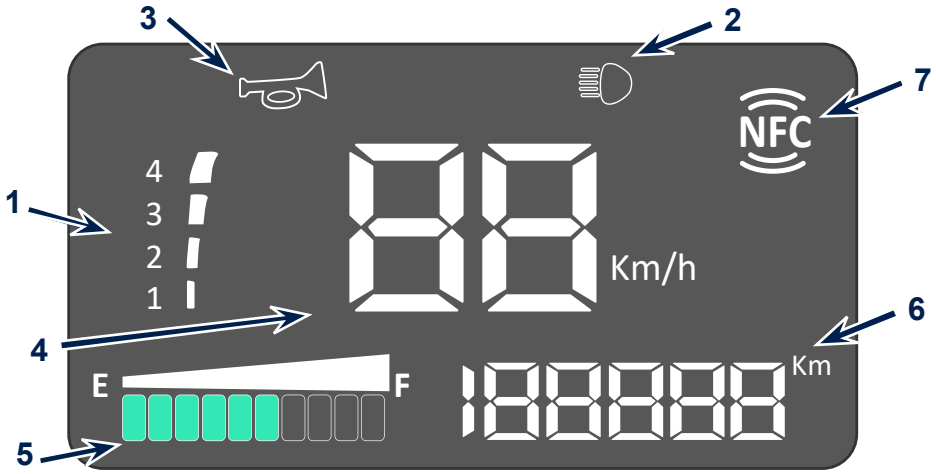
EN | The brake lever allows you to switch between electric and manual driving modes.
■ **Electric driving (D):** When the brake lever is in the rear, in the electric driving position, the scooter remains in a braking state and cannot be pushed. At this time, the scooter can be driven electrically in a normal way.
■ **Manual driving (N):** When the brake lever is forward, in the manual driving position, the scooter is in brake release state and can be pushed. This function can only be activated when the scooter is on and off. During this time, you cannot drive the scooter.

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EN | Before driving, familiarise yourself with local pedestrian traffic regulations and practice. Make sure the brakes are working properly, all functions are normal and the battery is fully charged before using it on the road.

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EN | 1. Sit on the scooter and press the power switch to turn it on. Swipe the NFC card over the control panel display to enable the scooter.
2. Adjust the speed control (4 levels) according to the road surface, surrounding conditions, and personal preference.
3. Hold the handlebar with both hands and slowly pull the throttle lever toward you with the fingers of your right hand to move forward.
4. When reversing, pay attention to your surroundings. Use the fingers of your left hand to slowly pull the throttle lever toward you.
5. To stop the scooter, gently release the throttle lever so it returns to the center position and the scooter stops.
6. To turn off the scooter, press and hold the power switch for 3 seconds, or swipe the NFC card over the display.

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EN | The factory NFC card is already paired and does not require pairing again. Up to 10 cards/phones can be paired at the same time. If you purchase a new NFC card or want to pair an NFC-enabled phone, follow these steps:
1. With the scooter powered on and unlocked, press and hold the speed control and the horn simultaneously for 5 seconds until you hear a beep.
2. Release the buttons. The NFC icon will flash on the control panel display. Hold the card/phone near the NFC reading area.
3. A new beep confirms pairing and the system returns to the main screen.

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EN | 1. Turn off the scooter by pressing the power switch, then remove the removable battery by pulling its handle.
2. Connect the charger to the battery charging port and plug it into a power outlet to start charging.
3. The battery level is shown by the battery's indicator light: red = charging required, green = fully charged.
4. The battery can also be charged without removing it from the scooter. While charging, the scooter will be disabled for safety to prevent movement.
5. When charging is complete, unplug the charger, press the power switch, and swipe the NFC card to start the scooter.

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EN | The scooter has an automatic sleep mode: if there is no operation for 20 minutes, the controller enters sleep mode. To continue using the scooter, turn it off and turn it on again.
The battery includes reverse polarity protection. It also has a 72-hour sleep mode: if it has not been used for that period, press the battery power switch to activate it and use the scooter normally.

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EN

KMINA - LIGHTWEIGHT SCOOTER BLUE	K10094
KMINA - LIGHTWEIGHT SCOOTER GRAY	K10095

MATERIALS
Magnesium and carbon fiber.

WARNINGS
■ Read and understand the entire manual before first use. Failure to follow warnings may cause personal injury.
■ Maintain a high level of attention while driving and comply with local pedestrian and traffic rules where applicable.
■ Use safety elements outdoors when necessary (reflective items, lights, safety symbols, etc.).
■ Do not drive under the influence of alcohol. Before driving, familiarize yourself with scooter operation and practice at low speed.
■ Consult a healthcare professional if you have health conditions, physical limitations, or take medication that may affect safe driving.
■ Ensure moving parts do not cause injury, especially if children or animals are nearby.
■ Do not tie straps to the scooter.
■ Do not place fingers between parts. Wearing gloves is recommended when handling scooter parts (extra caution with metal/heavy parts).
■ Do not expose the scooter to direct sunlight for prolonged periods, there is a risk of overheating surfaces and metal parts.
■ The scooter is designed for one passenger only, do not carry companions.
■ Do not use the scooter for purposes other than its intended use.
■ Do not tow the scooter. Towing may exceed the speed threshold and damage critical components.
■ Do not hang objects on the handlebar.
■ Do not smoke while seated on the scooter. Keep ashtrays away from the seat and fully extinguish cigarettes before disposal.
■ Sit as far back as possible on the seat when getting on/off. Sitting too far forward may cause falls.
■ Do not pull yourself up using the handlebar. Do not lean on the handlebar when standing up.
■ Do not use the armrests to support body weight.
■ Do not put weight on the footrest/foot area, and do not concentrate weight on one side.
■ Avoid body movements that shift the center of gravity while seated.
■ Avoid sudden turns. Reduce speed and maintain stability during turns and maneuvers.
■ Do not drive with a depleted battery (risk of being stranded). Use extreme caution on curbs, ramps, and inclines.
■ Do not approach curbs/ramps diagonally. Approach straight on so both front wheels contact at the same time.
■ When descending a ramp, if speed is faster than expected, release the throttle lever to stop, then resume descending smoothly.
■ Reverse only on flat surfaces. Do not reverse on slopes.
■ To prevent tipping and overload, do not carry heavy items exceeding a total of 6.8 kg (15 lb).
■ A small drop-off can block the front wheel and cause a sudden stop or tip-over.
■ Do not reach or lean out of the seat, this may damage the backrest and cause a fall.
■ Do not attempt to cross obstacles higher than 25 mm, risk of tip-over/injury.
■ Before starting to move, set the speed selector to the lowest level.
■ In an emergency stop, release the throttle lever and allow the scooter to come to a complete stop.
■ Avoid sudden braking or acceleration on slopes.
■ Do not attempt to cross obstacles while going up or down a slope. Approach in a straight line and do not stop halfway over the obstacle.
■ Driving in narrow spaces must be done at minimum speed and with caution.

■ Seat height: the higher the seat, the higher the risk of tipping (due to center of gravity changes).
■ Do not park the scooter on a slope.
■ Do not switch off the scooter while moving, unless there is no other option.
■ Do not get off the scooter on a slope or incline.
■ Do not drive on potentially dangerous roads or slopes: snow, ice, grass, leaves, or other slippery/unstable conditions.
■ Do not drive on uneven terrain or soft surfaces. Avoid tall grass that may become entangled in the mechanisms.
■ Avoid loose gravel and/or sand.
■ Do not use on or near railroad tracks and crossings.
■ Do not ride near edges of rivers, lakes, or the sea. Do not use to cross water.
■ Do not expose the scooter to open flames.
■ Do not take the scooter into a shower/bathtub/pool/sauna. Rain, snow, salt, sprays, humidity, and icy surfaces can damage components and rust the chassis.
■ Avoid prolonged exposure to extreme heat or cold, it may affect components and the battery.
■ Keep hands, clothing, and objects away from wheels and tires while driving.
■ Do not drag objects behind the scooter. Loose clothing or objects may get caught in wheels/moving parts.
■ Do not use the scooter in manual driving mode without an attendant.
■ Do not remain seated on the scooter when it is in manual driving mode.
■ Do not use manual driving mode on a slope, there is a risk of unintended movement.
■ Do not shift the lever into manual mode while driving.
■ Incomplete folding/unfolding may cause damage and/or injury. Keep a safe distance while folding/unfolding.
■ Avoid folding the scooter on steep slopes.
■ To carry the scooter up/down stairs: use extreme caution, fold the scooter, and use several people if necessary. Do not use escalators.
■ Do not use the scooter inside a moving vehicle.
■ When transporting the scooter, do not place it on a front seat. Secure the scooter so it cannot roll or move.
■ During transport, secure the scooter and battery. The manual indicates removing the battery and transporting it packaged separately, with no contact with flammable materials.
■ When the scooter is not in use, removing the battery is recommended.
■ Switch off the scooter before getting on or off.
■ If using an elevator, switch off the scooter and remain still to prevent unintended movement.
■ Do not exceed the maximum permitted user weight.
■ Check electrical components (corrosion, wear, damage), cables, and terminals, replace if necessary.
■ If any seating system parts are loose, or upholstery is worn/damaged, repair or replace immediately.
■ Do not modify seat specifications or add unapproved cushions/padding, this may create an unstable posture and cause falls.
■ Adding accessories may change weight, dimensions, and center of gravity. Use extra caution after any change.
■ Transmitting equipment may affect scooter control and cause unintended movement or brake release.
■ Do not turn on portable transmitters while the scooter is powered on.
■ If unintended movement or brake release occurs, switch off the scooter as soon as it is safe.
■ Avoid using the scooter near powerful transmitters (radio/TV). If very close to broadcast antennas, verify normal operation before use.
■ Avoid placing the scooter adjacent to or stacked with other electrical equipment. If unavoidable, verify correct operation. Keep portable RF equipment at least 30 cm away from the scooter and its cables.
■ Before each use, perform a basic check: tire condition, electrical connections (no corrosion, secure), connectors properly fastened, battery charged, brakes operational, and correct driving mode selected.
■ Suitable for air transport.

BATTERY CHARGING AND MAINTENANCE

■ Use only the supplied charger. Do not use other types of chargers.
■ Switch off the scooter before charging. Do not use the scooter while charging.
■ Do not keep the lithium battery at a low charge for long periods, it reduces battery life and may damage it.
■ Disconnect the charger when the battery is fully charged. Do not extend charging unnecessarily.
■ Keep the charger connector clean and dry. Avoid humidity.
■ Do not handle or repair a faulty battery. Authorized technical service must intervene.
■ Replace the battery only with a specified/authorized battery to avoid incompatibilities.
■ Do not connect external devices to the electronic system or use the scooter battery to power other devices.
■ Protect the charger from moisture and always charge in a dry environment.
■ Do not use the charger if it has been dropped or is damaged.
■ Do not touch the power outlet with wet hands.
■ Do not dispose of the battery in household trash or fire. Take it to a recycling center and cover terminals before disposal.

MAINTENANCE AND CLEANING

■ Check wheel wear regularly to see if replacement is needed.
■ Check that the seat is not loose or excessively wobbling before driving.
■ Check that the brake meets braking requirements. The scooter can be placed on an appropriate incline for further inspection.
■ Inspect the motor to see if replacement is needed, especially after one year of use.
■ Clean the scooter surface frequently. Do not wash the scooter directly with water to avoid permanent damage.
■ Clean the scooter with a damp cloth.
■ Develop good driving habits and do not operate the scooter beyond the permitted range.
■ Store the scooter in a dry, ventilated place. Do not place it in water or a very humid environment (humidity must be below 85%).
■ If the scooter has not been used for a long time, take dust-prevention measures. Keep the battery fully charged and charge it at least once a month.

PERFORMANCE PARAMETERS

■ Maximum speed: 6 km/h.
■ Maximum slope: 9°.
■ Braking method: Electromagnetic.
■ Lithium battery: 24V/11.6Ah.
■ Motor: 24V/250W.
■ Input charger: 100-240V/50-60Hz
■ Output charger: 24V/2A
■ Control panel: 24V/120A.

WARRANTY

IMPROVING MOBILITY DEVELOPMENT S.L. warrants that this device is free from defects.
■ The warranty for the motor, battery, charger, and control panel is valid for 1 year from the purchase date.
■ The product warranty remains valid for 3 years from the consumer's original purchase date.
The warranty does not extend to parts subject to natural wear that require periodic replacement.
■ This warranty does not cover:
Spare parts supplied by anyone other than IMPROVING MOBILITY DEVELOPMENT S.L. or the authorized distributor.
Any failure during the warranty period that is not caused by a defect in materials or workmanship, or if the failure is caused by misuse.
IMPROVING MOBILITY DEVELOPMENT S.L. is not responsible for self-damage, third-party damage, or property damage, whether accidental or due to bad faith.
IMPROVING MOBILITY DEVELOPMENT S.L. is not responsible for theft or loss of parts or the scooter.
Any modification carried out by the user or by third parties without the approval of IMPROVING MOBILITY DEVELOPMENT S.L. will void the warranty coverage.

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K10095

EN | POSSIBLE FAILURES AND SOLUTIONS

	FAILURES	POSIBLE CAUSES	SOLUTIONS
1	Control lever failure or low maximum Speed	Low battery voltage	Fully charge the battery
		Electromagnetic handle not engaged	Engage the electromagnetic brake
2	Motor does not work when power is turned on	Low battery voltage	Fully charge the battery
		Electromagnetic handle not engaged	Engage the electromagnetic brake
		Loose or faulty motor wire connector	Seek professional help or repair service
3	Insufficient driving range after a single charge	Inadequate charging or charger mal-function	Fully charge or replace the charger
		Excessive uphill driving, strong headwinds, frequent braking/starting, or overloading	Avoid such conditions
		Battery aging or damage	Replace the battery
4	Charger does not charge	Loose charger socket or poor plug-socket connection	Tighten loose screws and connections
		Blown fuse in the charger/e-bike charging plug	Replace the fuse
		Disconnected battery wiring	Reconnect the wiring
5	Unusual noise while driving	Loose screws/nuts	Tighten loose screws and nuts
		Worn wheel bearings	Replace the bearings
6	Other faults	Faults that cannot be resolved independently or diagnosed	Contact the supplier or repair service im-mediately. Do not disassemble components without authorization (warranty void if tampered with)

EN | MAINTENANCE CHECK

INSPECTION ITEMS	AT ANY TIME	WEEKLY	MONTHLY	SEMI-ANNUAL
Each part			X	
Turning,driving,-setting,dismantle-ment,etc		X		
Brake	X			
Connecting lines		X		
Battery condition	X			
Tire Condition			X	
Motor				X
Devices		X		
Cleanliness	X			

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